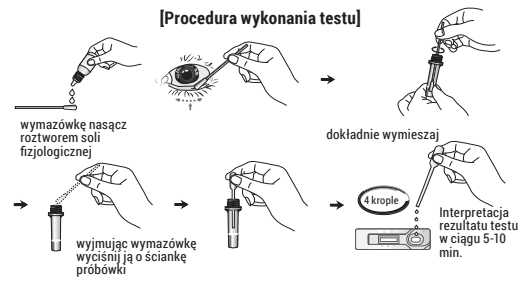
VetExpert Rapid Test CDV Ag

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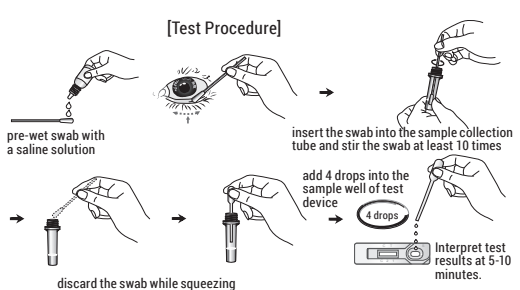
VetExpert Rapid CDV Ag Test

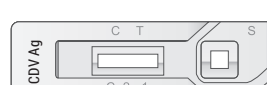
- Principles**
The VetExpert Rapid CDV Ag Test is a chromatographic immunoassay for the qualitative detection of Canine Distemper virus antigen in conjunctiva, urine, serum or plasma. The VetExpert Rapid Canine Distemper virus Ag Test has a letter of “T” and “C” as test line and control line on the surface of the device. Both the test line and control line in result window are not visible before applying any samples. The control line is used for procedural control. Control line should be always appeared. If the test procedure is performed properly and the test reagents of control line are working. A purple test line will be visible in the result window if there is enough Canine Distemper virus antigen in the specimen. The specially selected Canine Distemper virus antibodies are used in test band as both capture and detector materials. These enable the VetExpert Rapid CDV Ag Test to identify Canine Distemper virus antigen in conjunctiva, urine, serum or plasma with a high degree of accuracy.
- Contents of the test** (2, 5 or 10 tests depending on the kit)
VetExpert Rapid CDV Ag Test contains following items to perform the assay
 - 2,5 or 10 VetExpert Rapid CDV Ag Test
 - 2,5 or 10 tubes with thinner
 - 2,5 or 10 sticks for swabs
 - 2,5 or 10 disposable pipettes
 - Instruction for use
- Precautions**
 - For veterinary diagnostic use only.
 - For best results, strict adherence to these instructions is required.
 - All specimens should be handled as being potentially infectious.
 - Do not open or remove test kit from their individually sealed pouches until immediately before their use.
 - Do not use the test kit if the pouch is damaged or the seal is broken.
 - Do not reuse test kit.
 - All reagents must be at room temperature before running the assay.
 - Do not use reagents beyond the stated expiration date marked on the label.
 - The components in this kit have been quality control tested as standard batch unit. Do not mix components from different lot numbers.
- Storage and Stability**
The kit can be stored at room temperature (2-30°C) or refrigerated. The test kit is stable through the expiration date marked on the package label. DO NOT FREEZE. Do not store the test kit in direct sunlight.
- Specimen Collection and Preparation**
 - The test should be performed the conjunctiva, urine, serum or plasma.
 - After collecting the specimen using swab, the specimen should be immediately extracted and tested.
 - If specimens are not immediately tested, they should be refrigerated at 2-8°C. For storage not less than 48 hours, freeze the specimen at -20°C or below.

- Procedure of the test**
 - Collect the samples from conjunctiva or urine using the sample collection swab prewetted with saline solution (not provided). In case of serum or plasma samples, you can use the dropper.
 - Insert the swab into the specimen tube containing 300 µl of assay diluent. In case of serum or plasma samples, add 2-3 drops of the serum or plasma into the specimen tube containing 300 µl of assay diluent using the dropper.
 - Mix the swab samples with assay diluent to extract well.
 - Remove the test device from the foil pouch, and place it on a flat and dry surface.
 - Add four (4) drops of the mixed sample into the sample hole using the dropper, drop by drop and slowly
 - As the test begins to work, you will see purple color move across the result window in the center of the test device. If the migration has not appeared after 1 minute, add one more drop of the mixed sample to the sample well.
 - Interpret test results at 5-10 minutes.

[Test Procedure]

- Interpretation of the test**
A color band will appear in the left section of the result window to show that the test is working properly. This band is the control band. The right section of the result window indicates the test results. If another color band appears in the right section of the result window. This band is the test band.

- Negative result**
The presence of only one band within the result window indicates a negative result.

- Positive result**
The presence of two color bands (“T” and “C”) within the result window, no matter which band appears first indicates a positive result.

- Invalid result**
If the purple color band is not visible within the result window after performing the test, the result is considered invalid. The directions may not have been followed correctly or the test may have deteriorated. It is recommended that the specimen be re-tested.

- Limitation of the test**
Although the VetExpert Rapid Canine Distemper virus Ag Test is very accurate in detecting Canine Distemper virus antigen, a low incidence of false results can occur. Other clinically available tests are required if questionable results are obtained. As with all diagnostic tests, a definitive clinical diagnosis should not be based on the results of a single test, but should only be made by the veterinarian after all clinical and laboratory findings have been evaluated.

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